
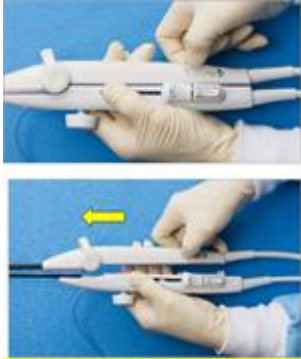

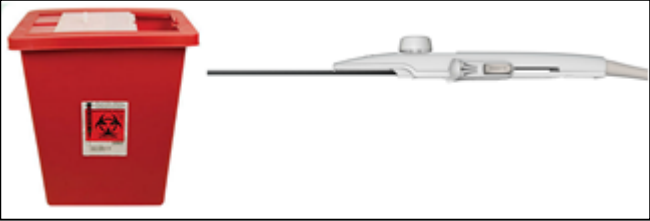




Post-Treatment Processes

Post-Procedure processes for the Support Staff include the following:

PROCESS	DESCRIPTION
Disconnect Dispersive Electrodes	<ul style="list-style-type: none"> • Disconnect the Dispersive Electrodes from the RF Generator. • Remove the Dispersive Electrodes from the patient. • Dispose of the Dispersive Electrodes. • Inspect the patient's skin for redness. 
Disconnect Sonata Treatment Device Components	<p>To separate the Treatment Device components:</p> <ol style="list-style-type: none"> 1. Hold the device with both hands. 2. Lift the Release Latches simultaneously.

PROCESS	DESCRIPTION
	<p>3. Slide the RFA Handpiece forward.</p> <p>4. Separate the IUUS Probe from the RFA Handpiece.</p> 
<p>Pre-Clean IUUS Probe</p>	<p>Clean off the visible blood and tissue with a towel or sponge moistened with water or enzymatic spray/solution. Then, place in the transport tray and send for Reprocessing.</p> <p>Refer to the <i>Sonata Operator's Manual</i> for complete instructions.</p> 

PROCESS	DESCRIPTION
<p>Dispose of RFA Handpiece</p>	<p>The RFA Handpiece is a single use device. Use caution when handling; treat as sharps. Dispose of the RFA Handpiece as a medical biohazard with sharps.</p> <p>Refer to the <i>Sonata Operator's Manual</i> or the <i>Post Procedure</i> section in the <i>Quick Reference Guide</i> for complete instructions.</p> 
<p>Video Recording Management</p>	<p>A video of the procedure is recorded automatically. A window displays when the IUUS Probe is disconnected in which the choice is given to Keep or Discard the recording:</p> <ul style="list-style-type: none"> • There is no SAVE button. • The Keep option saves the recording on the Sonata® System's hard drive. • The Discard option deletes the recording permanently. 

PROCESS	DESCRIPTION
Power OFF and System Care	<ul style="list-style-type: none"> • Power OFF the Ultrasound Console with the Power icon on the Ultrasound Display. • Power OFF the RF Generator by using the switch in the back. • Wipe the system, including the mouse and footswitch, with the low alcohol, medical-grade disinfectant. • Wrap the power cables for storage. 
Support Staff Responsibilities	<p>Support Staff responsibilities include, but are not limited to, the following:</p> <p>Setup</p> <ul style="list-style-type: none"> • Sterile the draped table. • Set up the Sonata[®] System by powering on the system and connecting all the cables. • Prepare the Patient by applying the Dispersive Electrodes.

PROCESS	DESCRIPTION
	<p>Procedure</p> <ul style="list-style-type: none"> • Click the Tip. • Freeze and store the image. • Press Ready after the Final Safety Check. <p>Post-Procedure</p> <ul style="list-style-type: none"> • Disconnect the cables. • Remove the Dispersive Electrodes from patient and inspect the skin. • Pre-clean the IUUS Probe for reprocessing. • Dispose of the RFA Handpiece. • Power OFF the Sonata® System. • Wipe down the system, including the Footswitch and Mouse. • Wrap and store the power cord. • Store the Footswitch. • Lower the System Cart completely.